

INSTRUMENT TABLE "ANNA"

OPERATING THEATRE



MADE IN SWEDEN

INSTRUMENT TABLE "ANNA"

- Flexible stainless steel tables of high quality in several unique designs
- Sound absorbing table top, easy to clean
- Electrical maneuvered with one single button, large height span
- Steady also on uneven floors

The instrument table called "Anna" are of highest quality and are available in several standard designs but can also be tailor-made. The table tops are sound absorbing and are provided with drip pan. The wheelbase is manufactured in stainless steel for easy cleaning.

Anna Tripod has a flat lifting column assemble at the short side of the table which allows great accessibility over operating table. The wheelbase with three wheels makes the table steady even on uneven floors and the location of the large wheels and their brake function makes the table extra steady. The height control can be placed where it gives optimal accessibility.

Anna Cross was developed in co-operation with operating personnel with need of a larger instrument table with high lifting capacity and good accessibility. The wheelbase consists of round, stainless steel pipes with four easily rolled, ball bearing and lockable wheels. The design of the base gives plenty of space for feet and chairs. The table-top has also a stainless steel underneath for easy cleaning.

All of Rini instrument tables fulfil the highly set demands in challenging operating theatres.



TECHNICAL INFORMATION

GENERAL

	<i>"Tripod"</i>	<i>"Cross"</i>
Length	730 mm	970 mm
Width	550 mm	620 mm
Lifting capacity	50 kg/110 lbs	75 kg/165 lbs
Weight unit	23 kg/51 lbs	43 kg/95 lbs
Height adjustment from floor	700-1100 mm	850-1350 mm

TABLE-TOP


	<i>"Tripod"</i>	<i>"Cross"</i>
Length	700 mm	1200 mm
Width	500 mm	600 mm

ELECTRICAL SPECIFICATION

Battery	24 V
Charger	100-240 V
Charge time	8-10 hours

Other

Table of Anna Tripod type is at inquiry available with four wheels.

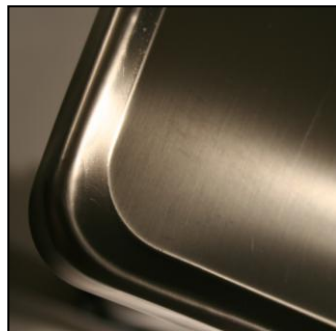

Rini Ergoteknik AB Fågelsångsvägen 3B
SE-186 42 Vallentuna Sweden
Tel +46 8-594 771 70 Fax +46 8-929 279
info@rini.se www.rini.se